8605 Santa Monica Blvd #48523 Los Angeles, CA 90069-4109

USA Phone: +1 (719) 203-1544 EU - Denmark Phone: 0045 8987 0700 info@endoca.com www.endoca.com ISO 14001: 2004 certified; ISO 9001: 2008 certified HACCP certified; GMP certified

Responsible Supervisor:	Martin V.
Sample _	Batch# 589
Date samples received:	13-February 2017
Date analysis began: 🔤 🗌	13-February 2017
Date sample report produced: _ ID Number when available: _	14-February 2017
Sample Mass	1 g


Endoca 1.95% Total CBD: **Cannabinoid Profile**

HPLC Chromatograph Raw Data

Component	Mass (%)	Amount (mg/g)	Limit	0.00 2487Channel 1 - 225 nm 0.00 B B B B B B B B B B B B B B B B B B B
CBD	1.73	17.30	N/A	P 1 -0.20 ℃ ↓
CBDA	0.22	2.20	N/A	
CBG	0.02	0.20	N/A	2487Channel 2 - 300 nm
CBGA	<0.01	<0.10	N/A	0.002
CBN	0.07	0.70	N/A	₹ 0.00
D9THC	<0.20	<2.00	N/A	
THCA	<0.20	<2.00	N/A	2.00 4.00 6.00 8.00 10.00 12.00 14.00 16.00 18.00 Minutes
Total CBD	1.95	19.50	N/A	



Cannabinoids as Percent of Total Mass

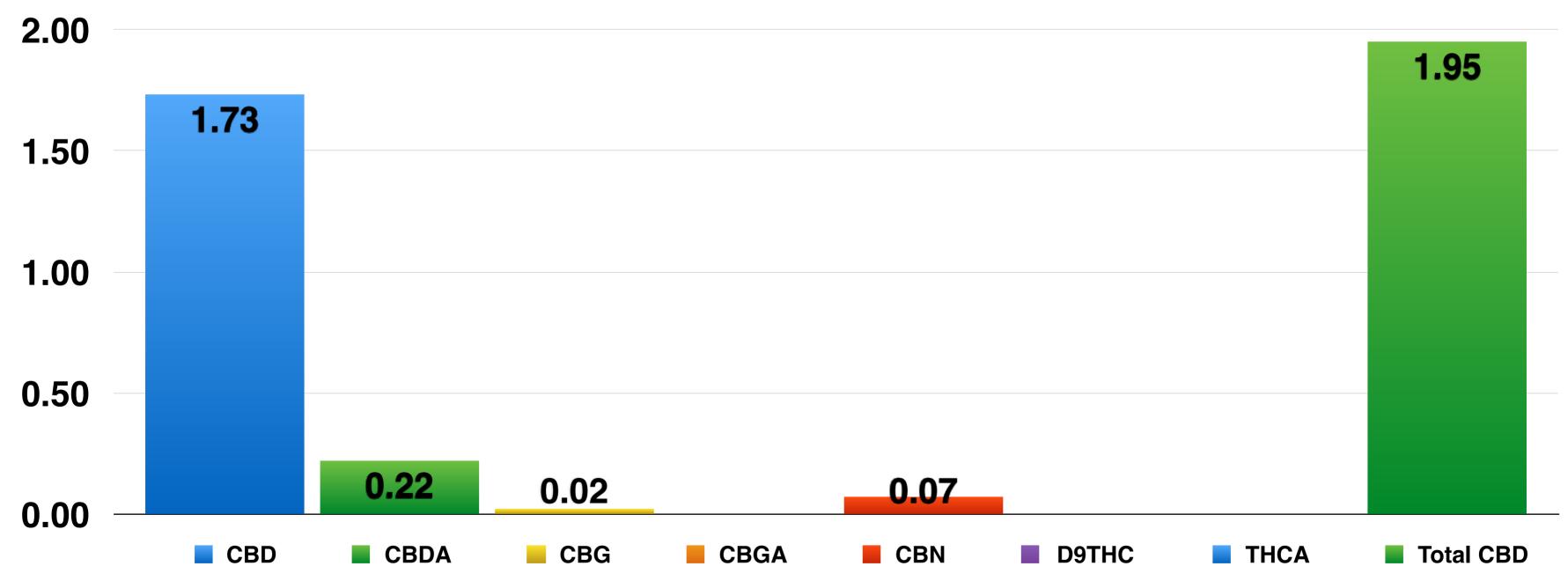

This report and all information herein has been confirmed by a certified lab and shall not be reproduced, except in its entirety, without the written consent of Endoca. The report is for informational purposes only and should not be used to diagnose, treat, or prevent any medical related symptoms. Results are applicable only for samples tested, and for the specific tests conducted.

All tests carried are done so under strict certified laboratory protocols, guidelines and supervision. See appendix for more information or contact us at <u>www.endoca.com</u>


8605 Santa Monica Blvd #48523 Los Angeles, CA 90069-4109


USA Phone: +1 (719) 203-1544 EU - Denmark Phone: 0045 8987 0700 info@endoca.com www.endoca.com ISO 14001: 2004 certified; ISO 9001: 2008 certified HACCP certified; GMP certified

Responsible Supervisor: Martin V. Sample **Batch# 589 Date samples received: 13-February 2017** Date analysis began: **13-February 2017** Date sample report produced: **14-February 2017 ID** Number when available: Sample Mass **1** g

Endoca 1.95% Total CBD: **Terpenoid Profile:**

Component	Amount %	FID1 A, (TERPENE 2017-02-14 11-46-16\002F0201.D)
β-Caryophyllene	0.29	
a-Humulene	<0.01	40 - 66 6
Caryophyllene oxide	<0.01	23.544
Myrcene	0.28	35 -
a-Pinene	0.57	21.083
Terpinolene	<0.01	30 -
Humulene epoxide II	<0.01	5
Limonene	0.77	25 - 25-
β-Pinene	<0.01	53 A410 (A10 (A10 (A10 (A10 (A10 (A10 (A10 (A
	.0.04	

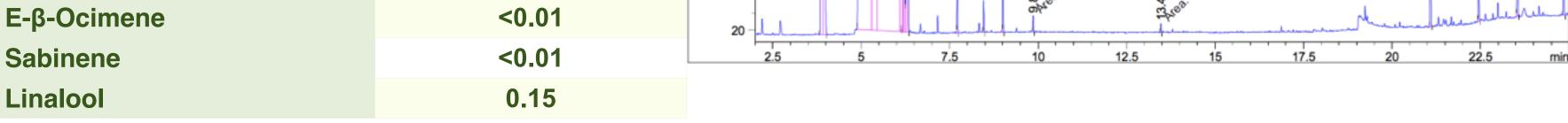

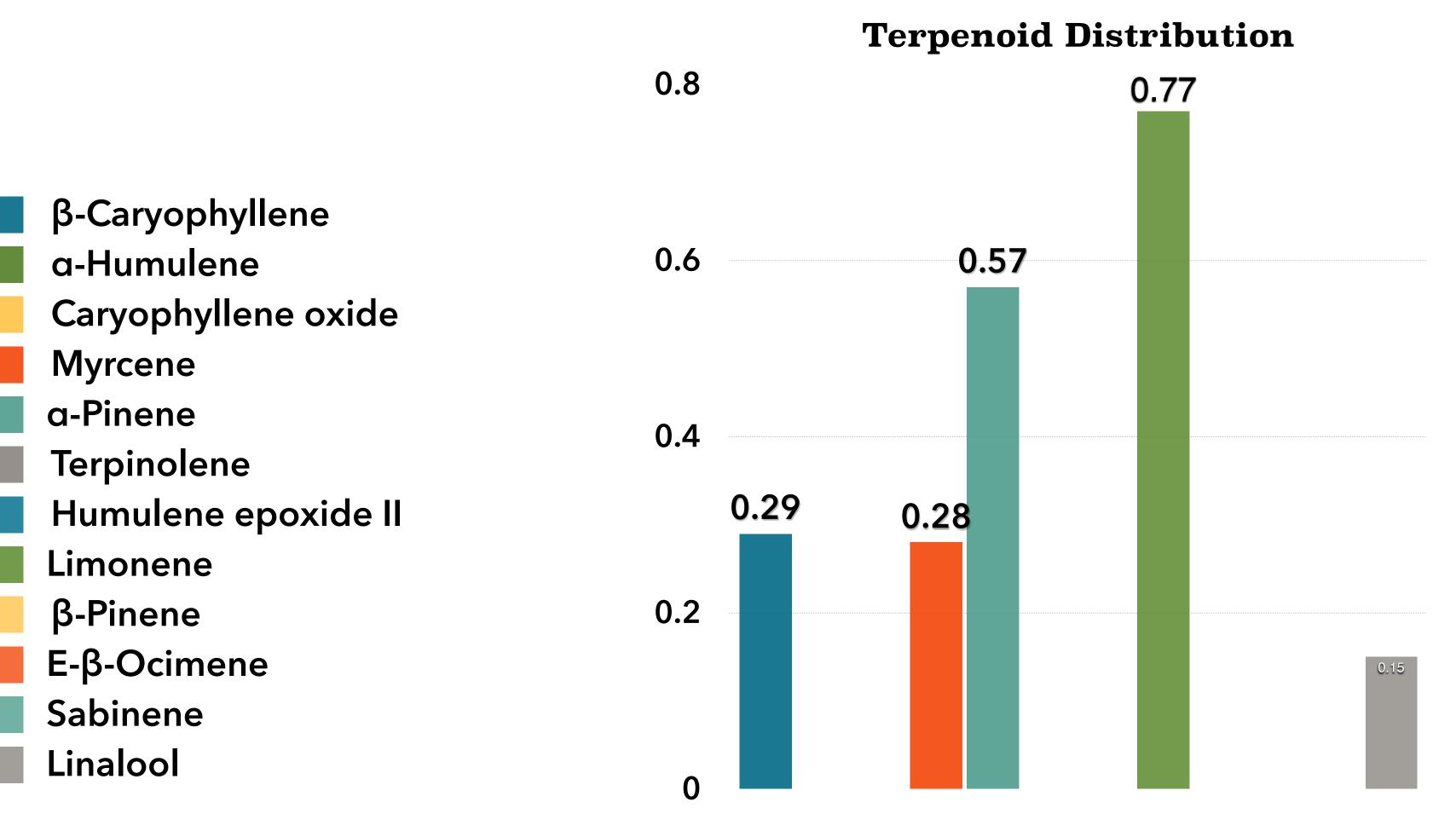

This report and all information herein has been confirmed by a certified lab and shall not be reproduced, except in its entirety, without the written consent of Endoca. The report is for informational purposes only and should not be used to diagnose, treat, or prevent any medical related symptoms. Results are applicable only for samples tested, and for the specific tests conducted.

All tests carried are done so under strict certified laboratory protocols, guidelines and supervision.

See appendix for more information or contact us at <u>www.endoca.com</u>

5

Amount (%) Mass of Mycotoxins

0

Α

B1

8605 Santa Monica Blvd #48523 Los Angeles, CA 90069-4109

USA Phone: +1 (719) 203-1544 EU - Denmark Phone: 0045 8987 0700 info@endoca.com www.endoca.com ISO 14001: 2004 certified; ISO 9001: 2008 certified HACCP certified; GMP certified

Responsible Supervisor:

Sample **Date samples received:** Date analysis began: Date sample report produced: ID Number when available: Sample Mass

Martin V.

1 g

Batch# 589 13-February 2017 13-February 2017 14-February 2017

Endoca 1.95% Total CBD: **Microbial Profile:**

Component		Results
Listeria m.	1 g	Not Detected
Escherichia c.	1 g	Not Detected
Salmonella	25 g	Not Detected
Yeast	1 g	<1 x 10 ²
Mould	1 g	<1 x 10 ²

All Mycotoxins at Non Detectable (ND) levels

NO MYCOTOXINS RESIDUES TO DISPLAY

B2

G1

G2

Nutrition Facts

Component	%
Moisture	<0.1
Protein	ND
Total fat	ND
Total Carbohydrates	ND
Dietary Fibers	ND
Sugars	ND
Ash	ND

This report and all information herein has been confirmed by a certified lab and shall not be reproduced, except in its entirety, without the written consent of Endoca. The report is for informational purposes only and should not be used to diagnose, treat, or prevent any medical related symptoms. Results are applicable only for samples tested, and for the specific tests conducted.

All tests carried are done so under strict certified laboratory protocols, guidelines and supervision. See appendix for more information or contact us at <u>www.endoca.com</u>



8605 Santa Monica Blvd #48523 Los Angeles, CA 90069-4109

USA Phone: +1 (719) 203-1544 EU - Denmark Phone: 0045 8987 0700 info@endoca.com www.endoca.com ISO 14001: 2004 certified; ISO 9001: 2008 certified HACCP certified; GMP certified

Responsible Supervisor:

Sample Date samples received: Date analysis began: Date sample report produced: **ID** Number when available: **Sample Mass**

Martin V. **Batch# 589**


<u>13-February 2017</u> 13-February 2017 14-February 2017

Endoca 1.95% Total CBD: **Heavy Metals Profile:**

All Heavy Metals at Non Detectable (ND) levels

1 g

Component	Mass (%)	Amount (mg/g)	Limit
Arsenic	< 0.01	ND	ND
Cadmium	< 0.01	ND	ND
Lead	< 0.01	ND	ND
Mercury	< 0.01	ND	ND
Flammables	< 0.01	ND	ND

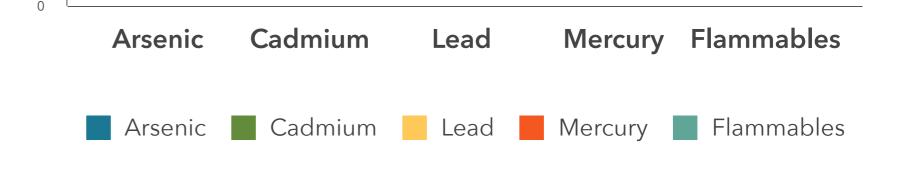

Conclusions:

No heavy metal residues detected.

No flammable residues detected.

No chemical residues detected.

This report and all information herein has been confirmed by a certified lab and shall not be reproduced, except in its entirety, without the written consent of Endoca. The report is for informational purposes only and should not be used to diagnose, treat, or prevent any medical related symptoms. Results are applicable only for samples tested, and for the specific tests conducted

All tests carried are done so under strict certified laboratory protocols, guidelines and supervision. See appendix for more

information or contact us at www.endoca.com

8605 Santa Monica Blvd #48523 Los Angeles, CA 90069-4109

USA Phone: +1 (719) 203-1544 EU - Denmark Phone: 0045 8987 0700 info@endoca.com www.endoca.com ISO 14001: 2004 certified; ISO 9001: 2008 certified HACCP certified; GMP certified **Responsible Supervisor:**

Sample Date samples received: Date analysis began: Date sample report produced: ID Number when available: Sample Mass Martin V.

Batch# 589 13-February 2017 13-February 2017 14-February 2017 1 g

Pesticide Analysis: Our tests looked for residue of

nearly 300 known pesticides finding no evidence

of any over detectable limits.

Endoca Labs tests our products thoroughly. Nearly 300 of the below pesticides concentrations were measured and we are proud to say that all tests measured below our detectable limits. Most tests have a threshold of 0.01 mg/k, while only a handful of tests have a threshold value of <0.05 mg/kg. Not a single test of Endoca products went over detectable threshold limits.

PESTICIDES MEASURED

Acrinathrin Azoxystrobin Biphenhin Bitertanol Biphenyl Bromopropylate Bromuconazole Bupirimate Cadusafos Captafol Captan Chlorphenson Chlorfenapyr Chlorfenvinphos Chlorothalonil Chlorprophame 3,5-Dichloraniline Chlorpyriphos Chlorpyrifos-methyl Chlorthal-dimethyl Cyfluthrin Cypermethrin Cyproconazole Cyprodinil Clomazone o,p-DDE P,P-DDE o,p-DDD P,P-DDD o,p-DDT p,p-DDT Deltamethri Diazinon Diclofop-methyl Dieldrin Dichlobenil Dichlofluanid Dichlorvos Dicloran Dicofol Dicrotophos Diethofencarb Diflubenzuron Dimetachlor Diniconazole Dodemorph Diphenylamine Alpha-Endosulfan Beta-Endosulfan Endosulfan-sulphate Ethion Etofumesate Ethoprophos Ehtoxyquin Etoxazole Etridiazole Etrimphos Famoxadone Fenarimol Fenazaquin Fenchlorphos Fenhexamid Fenihothion Fenpropidin Fenpropimorph Fenvalerate Formothion Fipronil Fipronil-sulfone Fludioxonil Flusilazole Flutriafol Folpet Fuberidazole Furathiocarb Hexaconazole HCB Alpha-HCH Beta-HCH Delta-HCH Heptachlor Heptachlor-epoxidceis Heptachlorepoxidtreans Iprodione Iprovalicarb Lambda- cyhalothrin Lindane Mecarbam Metalaxv Metazachlor Methidathion Metribuzin Mevinphos Myclobutanil Nuarimol Orthophenylphenol Oxadixyl Paclobutrazol Parathion Parathion-methyl Paraoxonmethyl Paraoxon-ethyl Penconazole Pendimethaline Permethrin Phenthoate Phorate Procymidone Profenofos Propiconazole Propyzamide Pyrazophos Pyrethrins Pyridaben Pyrimethanil Pyriproxyfen Quinoxifen Quitozene Pentachloraniline Phosphamidon Pyrifenox Prometryn Propanil Propoxur Proquinazid Prothiofos Simazine Spiroxamine T au-fluvalinate T ebuconazole T ebufenpyrad T ecnazene T efluthrin T erbuthylazine T etraconazole T etradifon T etramethrine T olclofos-methyl T olylfluanid Transfluthrin Triadimephon Triadimenol Trialate Trifloxystrobin Triflumizole Vinclozolin DDT isomersum Heptachlor (heptachloarnd heptachloer poxidsum) Trifluraline Chlorobenzilate 3-Chloraniline Abamectin (AvermectinBla and AvermectinBlb sum) Acetamiprid Aldicarb Aldikarbsulphone Aldicarbsulphoxide Azinphos-ethyl Azinphos-methyl Benalaxyl Benfuracarb Boscalid Buprofezin Carbaryl Carbendazim Carbofuran 3-hydroksicarbofuran Carbosulfan Chloridazon Cymoxanil Clofentezin Clothianidin Demeton-S-methyl Demeton-S-methyslulfoxid Diafenthiuron Difenoconazole Dimethoate Dimethomorph Diuron EPN Epoxiconazole Ethirimol Etofenprox Fenamidone Fenbuconazole Fenbutatinoxid Fenoxycarb Fenpyroximate Fenpropathrin Fensulfothion Fenthion Fenthionsulphone Fenthionsulphoxide Fluazinam Flufenoxuron Fluquinconazole Fonofos Formetanate Fosthiazate Hexythiazox Imazalil Imidacloprid Indoxacarb Isofenphos Methacrifos Isofenphosmethyl Krezoxim-methyl Linuron Lufenuron Malaoxon Malathion Mepanipirim Mepronil Metamitron Metconazole Methamidophos Methiocarb Methiocarbsulphone Methiocarbsulfoxide Methomyl Methoxyfenozide Metobromuron Monocrotophos Monolinuron Omethoate Oxamyl Pencycuron Phenmedipham Phosalone Phosmet Phosmeot xon Phoxim Pymetrozine Piperonylbutoxide Pyraclostrobin Pyridaphenthion Pyridate Pyrifenox Pirimicarb Pirimicarbdesmethyl Pirimiphos-methyl Primisulfuron-methyl Prochloraz Propamocarb Propargite Prothioconazole Prothioconazole-desthio Quinalphos SpinosynA SpinosynD Sulfotep T ebufenozide T eflubenzuron Thiabendazole Thiacloprid Thiamethoxam Thiodicar Thiophanatemethyl Tralkoxydim Triazophos Trichlorfon Triflumuron Triforine Triticonazole Zoxamide Acephate Amitraz Fenamiphos Fenamiphosulphone Fenamiphosulfoxid Nitempiram Fenthionoxonsulphone Fenthionoxonsulfoxid Kumapho Piriphenox Mehibuzine DEET

Our laboratory analysis is standardized after following protocols: LST EN ISO 6579:2003 / AC:2006 / P:2007 LST EN ISO 11290-1:2003 / A1:2004 / P:2005 LST ISO 16649-2:2002 / P:2009 LST ISO 21527-2:2008 Method PLM 486G

Note on Cannabinoid Testing:

All cannabinoids in their acid forms (ending in "-A") are convertible to their non-acid forms via a decarboxylation process (heating). The components lose mass through this process. To find the total theoretical active cannabinoids, one multiplies the acid forms by 87.7%. For example, CBD-A can be converted to active CBD using the formula: CBD-A x 0.877 = THC. In this

case, the Max CBD for the sample is: Max CBD = $(CBD-A \times 0.877) + CBD$. This method has been validated according to the principles of the International Conference on Harmonisation.

Chromatographic Analysis:

Analysis of cannabinoids content was performed using Waters 2695 (Milford, MA, USA) separation module equipped with auto injector, sample cooler, vacuum degasser and column heater units. Separation of all cannabinoids was accomplished on YMC PRO C18 (150 x 4 mm I.D., S-3 μ m) RP column coupled with C18 precolumn maintained at 30 °C by a CTO-20AC column oven.

Isocratic elution consisted of acetonitrile:water (4:1) was done in 30 min. The flow rate was maintained at 0.8 ml/min. The cannabinoids CBD, CBG and THC were monitored at 225 and CBDA, CBGA and THC-A were monitored at 300 nm respectively using dual absorbance detector Waters 2487 (Milford, MA, USA). The injection volume of 0.1 mg/ml sample was 10 μ l. Data evaluation was performed using Clarity software.

Quantification of cannabinoids was obtained from linear regression equation of calibration curve of individual reference standard by plotting concentration versus the area ratio.

Sample preparation for HPLC analysis

0.01 g (\pm .0001) of homogeneous cannabis extract was diluted with 1 ml of methanol (HPLC grade). Solution was sonicated for 5 min and vortexing for 10 sec. Samples before HPLC analysis were centrifuged at 4800 rpm and further diluted with methanol to the final concentration of 1 mg/ml.

Analysis of terpenes was performed using GC-FID system equipped with auto injector. Separation was accomplished on RTX-5 w/Integra-Guard, 30m, 0.25 mm ID column.

This report and all information herein has been confirmed by a certified lab and shall not be reproduced, except in its entirety, without the written consent of Endoca. The report is for informational purposes only and should not be used to diagnose, treat, or prevent any

medical related symptoms. Results are applicable only for samples tested, and for the specific tests conducted.

All tests carried are done so under strict certified laboratory protocols, guidelines and supervision. See appendix for more information or

contact us at www.endoca.com